


# A&D Pulse Ox Client Guide

## Instructions For Use


### 1 Insert Finger Into Pulse Ox

Insert **one finger** completely into the measurement section of the Oximeter with your **finger nail pointing up** towards the screen.

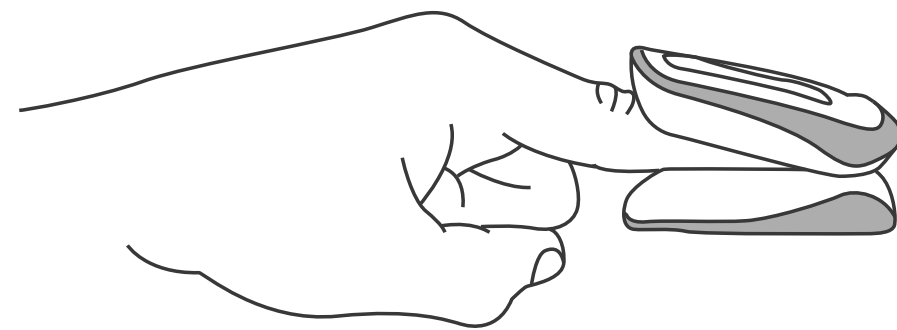


### 2 Press Power On Pulse OX

Stand within **20 feet** of the Tenovi Gateway and **press the power button** on the oximeter.

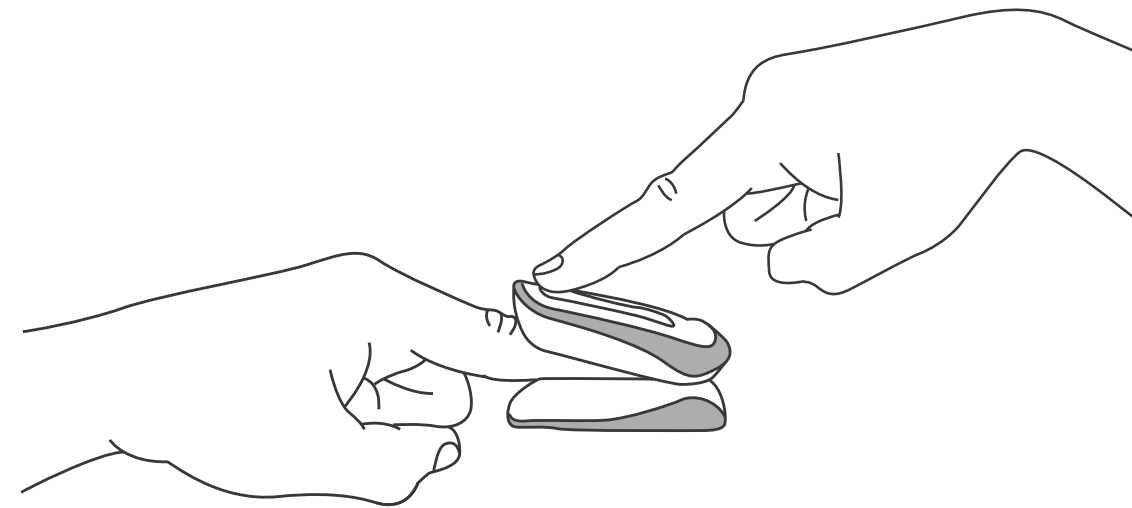


### 3 Wait 20 seconds

Keep your finger in the oximeter for **20 seconds** or until you hear a beep from the Tenovi Gateway. Then, **remove your finger** and wait for Tenovi Gateway to turn **GREEN**.

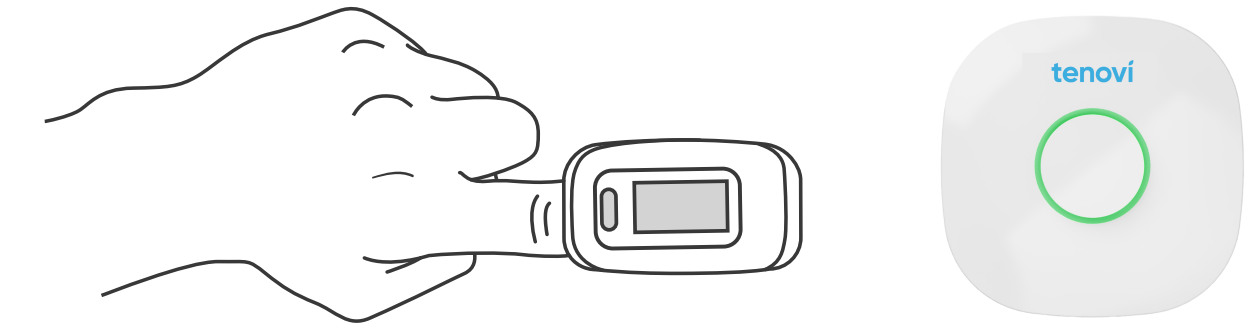

# A&D Pulse Ox Client Guide

## Device Diagram

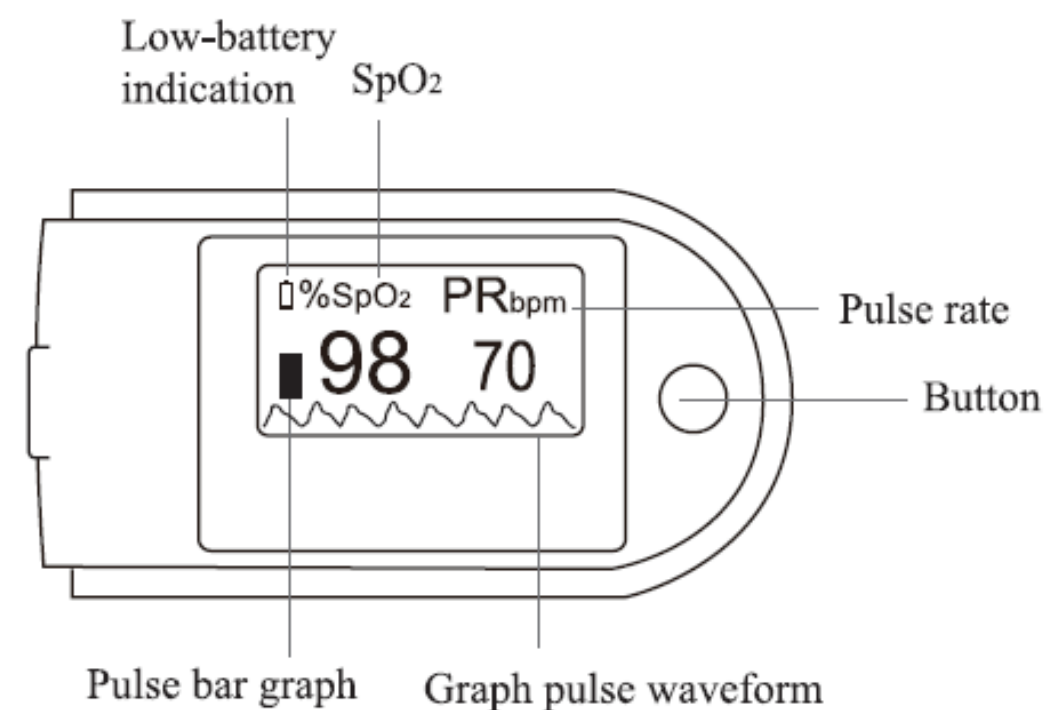

## Device Measurement Quirks and Notes

- Must remove finger from pulse ox for the gateway to turn green and transmit the measurement.
- If the user takes a second measurement before the device has turned off, the Pulse Ox will record the measurement in its storage but may fail to send it to the Gateway. However, that second reading will be sent to the gateway on the next connection cycle. We recommend waiting for the Pulse Ox to turn off before taking a second measurement.

- This oximeter is not intended for continuous monitoring.
- Measurement accuracy may be influenced by fingers that are too cold, too thin, or if fingernails are too long. You can use a larger finger such as your thumb or middle finger
- To avoid affecting its performance DO NOT place the oximeter in direct sunlight, high temperature, humidity, dust, cotton wool or near water. If the oximeter gets wet, please stop using it.
- When the oximeter is carried from cold environment to warm or humid environment, please wait 4 hours before using.

## Maintenance

- There is no need to calibrate the oximeter.
- Remove the batteries before cleaning.
- To clean, use 75% isopropanol alcohol dampened on a soft cloth. Then, wipe the oximeter for disinfection. Let air dry or dry with a soft cloth.
- DO NOT spray liquid directly on the oximeter