
1

Power Your Gateway

Plug your Gateway into an outlet with the provided power cord and wait until the LED on the Gateway turns red.



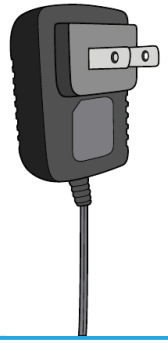

3

Use Your Blister Pack Daily

Your Smart Blister Pack will store open events. When you open a compartment while near the Gateway, the Gateway LED will turn green and all stored events will be sent to your care provider.



Measurement Needed



Measurement Successful

For further assistance, please visit Tenovi.com/support.

tenovi

© 2026 Tenovi Co.
All rights reserved. Republication, duplication or modification in whole or in part is prohibited without prior written permission from Tenovi, Co. All trademarks are owned by or licensed to Tenovi, Co. unless otherwise noted.
Revision 10.26.2026

2


Take Your Medication

Take your medication as advised by your doctor. Open the compartment corresponding to the current day and time. Push the pill through from the back or use the provided tool to puncture from the front.

Remove Your Tablet

Method 1

Push the tool through the compartment foil backing.

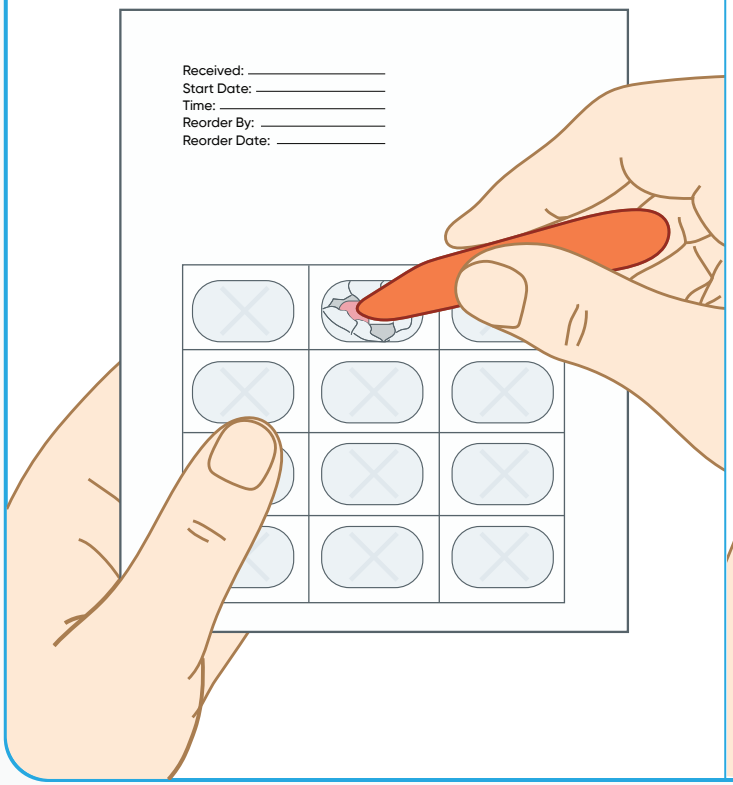

Method 2

Push the plastic bubble to force the medication through the foil backing.

